



MODfinder

Genetically Modified Organisms

Real-Time PCR detection and quantification



GENERON
Food Safety in our Genes

GMO is any living being added with an “extra-piece” of DNA in a lab

Using genetic engineering or transgenic technologies we can create combinations of plant, animal, bacterial and virus genes that can not occur in nature or through traditional crossbreeding methods

Real-time PCR is the gold standard in GMO analysis according to ISO norms which are now recognized worldwide as a benchmark

European Union was the first globally to introduce a legislation on the approval, traceability and detection of GMOs, including labelling of food and feed containing GMOs.

This provided maximum protection of public health and of environment, while at the same time providing a science based regulatory structure where biotechnology can flourish.

The development and application of reliable detection and quantitative analytical methods was essential for the implementation of the labelling rules.



Molecular analysis required to comply with the current European Union GMO legislation consist of three distinct steps:

- Detection
- Identification
- Quantification

Using MODIfinder products, along with an appropriate DNA extraction method, enables to test the presence of GMOs in grain commodities ingredients, semifinished and finished products, according to EU norms.

Insertion is based on a circular cloning vector with several elements

When the plasmid integrates into the host genome it takes a linear structure that accounts for the different PCR-based GMO tests to detect it These can be grouped into at least four categories corresponding to their level of specificity

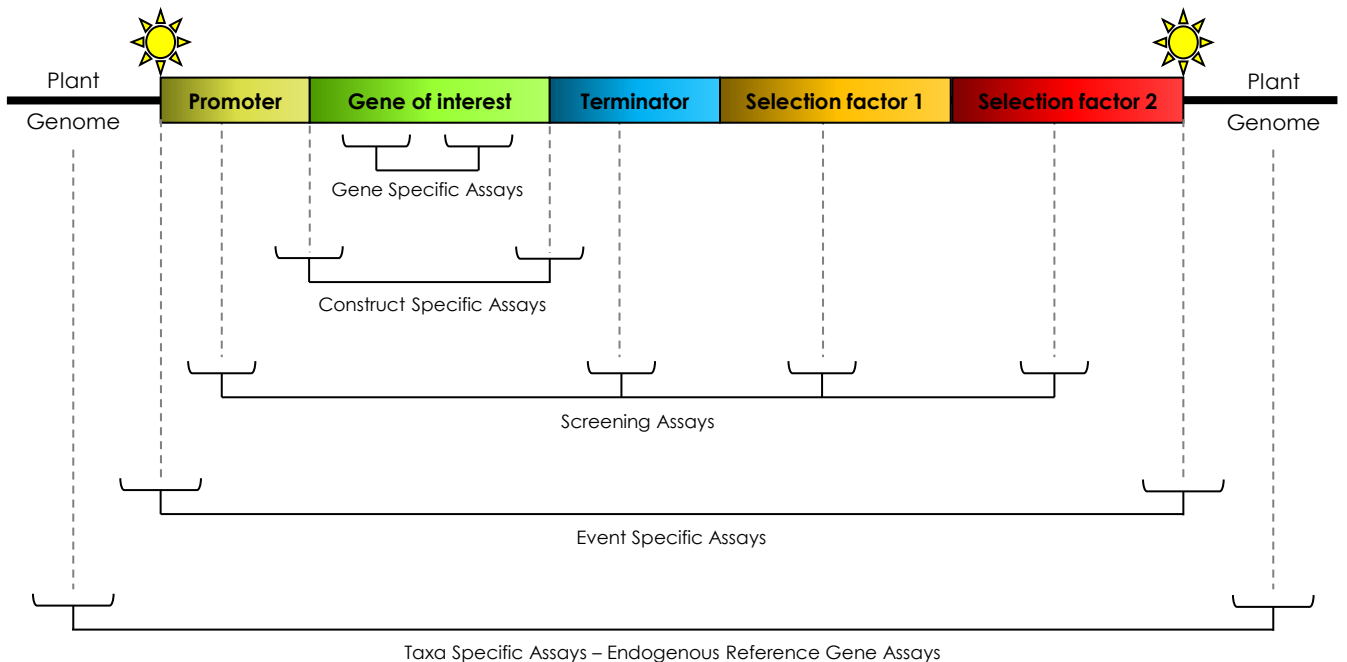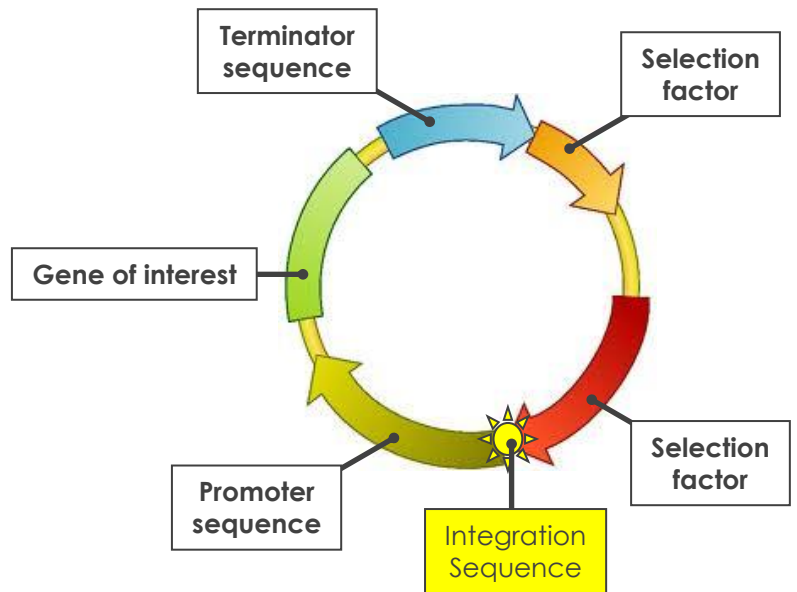

New GMO are constantly appearing and listed in databases


EU legislation covers all GMOs that have received Community authorization for the placing on the market; in contrast, unapproved GMOs are not permitted at any level.

Regulation (EC) 1829/2003 calls for the labelling of food and feed products that intentionally contain authorized GMO or have an adventitious (accidental and technically unavoidable) presence of authorized GMOs, above a threshold of 0.9% in otherwise non-GMO food and feed. When the sample contains a mixture of GMOs, each ingredient (analytically translated as content per taxon) may not exceed 0.9%, irrespective of its proportion in the final product.

As the number of genetically modified plants is increasing, it is necessary to use up-to-date methods for GMO screening and identification in food products

The basis for a qualitative GMO screening procedure is the use of control elements such as promoters, transcription terminators, and markers, such as resistance genes.

Negative responses from such a panel of screening methods eliminate the possibility of GMO presence in a test sample, but only if the selected screening methods cover all the GMOs to be detected.

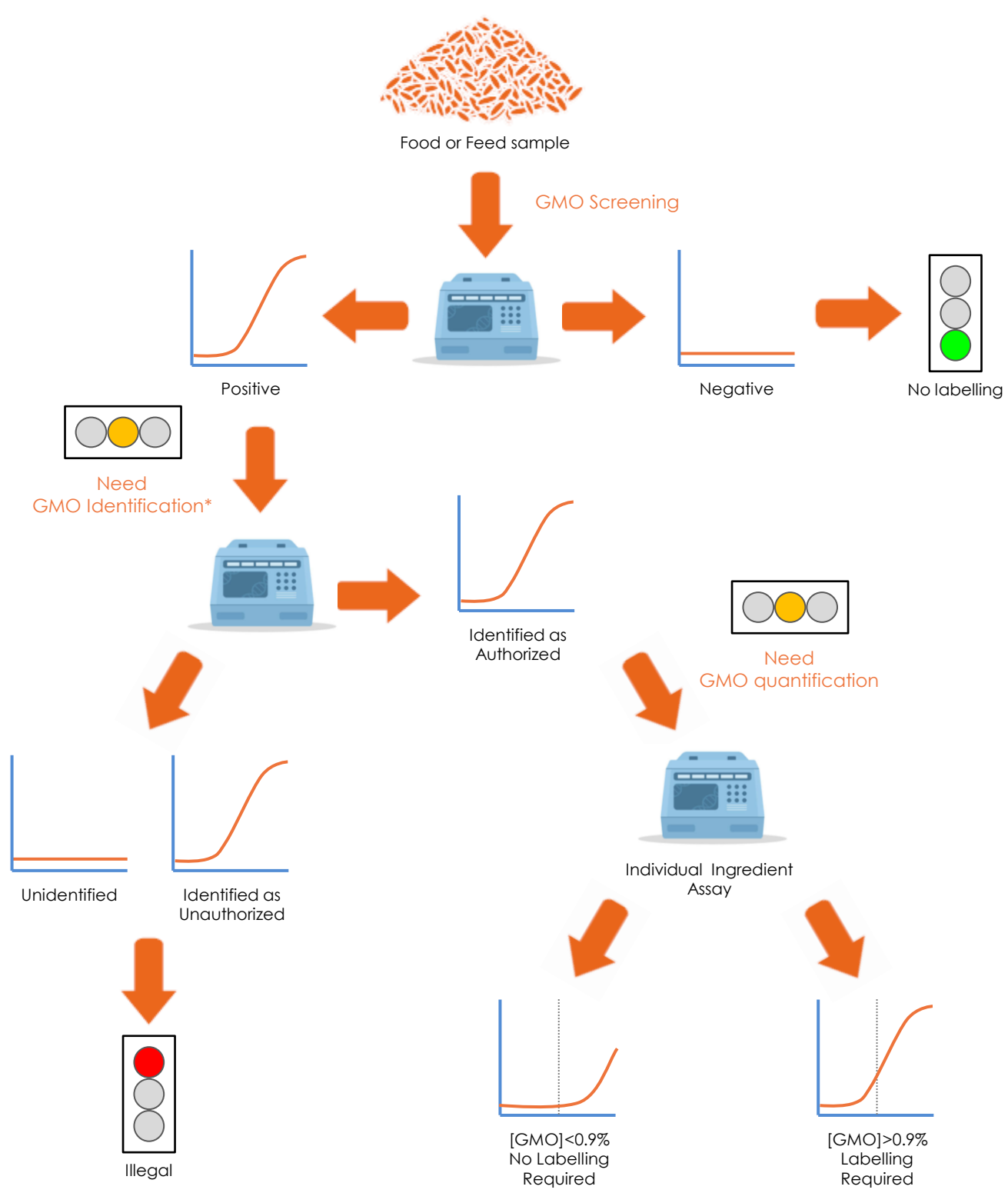

(*) When testing complex matrices composed by more then one ingredient is fundamental to have knowledge of the different plant varieties (taxa) present in the sample in order to restrict the number of GM specific test to be performed.

***MODfinder screening kits
detect the presence of
GMOs in a DNA extract***

Real-Time kits for all the relevant screening target (including gene and construct specific assays) allowing their detection in single-plex or multiplex. were developed following indications published by the EURL GMFF or in other peer reviewed scientific literature. Generon collection includes also Plant/Eukaryotic generic marker to determine the overall amplifiability of the extracted DNA according to ISO 24276.

Generon developed a portfolio of kits to provide the customers with a turn-key solution to test the presence of GMO according to EURL and ISO indications. The portfolio includes all the reagents and controls to execute the detection protocols



***MODfinder identification
kits detect the unique
signature of a
transformation event that is
the junction at the
integration locus between
the recipient genome and
the inserted DNA.***

The identification of the GMO contaminating the matrix under investigation is a fundamental step when considering the labelling rules enforced in Europe. Quantitation of GMO contamination must be in fact relative to something. In European legislation it refers to the presence w/w relative to the ingredient and not to the whole product. Therefore, knowing exactly which is the event to quantify is crucial. All the kits allow event detection with the sensitivity indicated by the method developers.

A. Singleplex Real-Time PCR kits for GM markers detection (50 reactions)

MODfinder 35S Promoter (p35S)	PGE01A-50	MODfinder 35S Terminator (t35S)	PGE03A-50
MODfinder NOS Terminator (tNOS)	PGE04A-50	MODfinder nptII	PGE07A-50
MODfinder PAT	PGE14A-50	MODfinder BAR	PGE15A-50
MODfinder pan-EPSPS	PGE16A-50	MODfinder FMV Promoter (pFMV)	PGE17A-50
MODfinder CryIAb/Ac	PGE25A-50	MODfinder NOS Promoter (pNOS)	PGE27A-50
MODfinder tE9 marker	PGE37A-50	MODfinder pSSuAra	PGE55A-50 (coming in 2023)

B. Multiplex* Real-Time PCR kits for GM markers detection (50 reactions)

MODfinder MultiSCREEN 2-plex p35S/tNOS	PGE05A-50
MODfinder MultiSCREEN 3-plex p35S/tNOS/pFMV	PGE26A-50
MODfinder MultiSCREEN 4-plex p35S/tNOS/pFMV + IAC	PGE26A-C-50
MODfinder MultiSCREEN 3-plex nptII/PAT/pan-EPSPS	PGE32A-50
MODfinder MultiSCREEN 4-plex p35S/tNOS/pFMV/BAR	PGE34A-50
MODfinder MultiSCREEN 4-plex nptII /PAT/pan-EPSPS /pNOS	PGE35A-50
MODfinder MultiSCREEN 2-plex nptII/pFMV	PGE40A-50
MODfinder MultiSCREEN 4-plex pan-EPSPS/PAT/BAR/nptII	PGE42A-50
MODfinder MultiSCREEN 2-plex pFMV/PAT	PGE46A-50
MODfinder MultiSCREEN 4-plex p35S/tNOS/pFMV/t35S	PGE47A-50
MODfinder MultiSCREEN 2-plex t35S/PAT+BAR	PGE51A-50
MODfinder MultiSCREEN 2-plex t35S/pFMV	PGE52A-50
MODfinder MultiSCREEN 2-plex pan-EPSPS/PAT	PGE53A-50
MODfinder MultiSCREEN 2-plex nptII/BAR	PGE54A-50
MODfinder MultiSCREEN 4-plex p35S/tNOS/PAT/BAR	PGE56A-50 (coming in 2023)
MODfinder MultiSCREEN 3-plex CryI Ab\Ac/tE9/mEPSPS	PGE57A-50 (coming in 2023)

C. Kits for DNA amplifiability verification (50 reactions)

MODfinder 28S (eukaryotic DNA)	PGE06A-50	MODfinder Chloroplast	PGE08A-50
MODfinder Soy	PGE09A-50	MODfinder Flax	PGE12A-50
MODfinder Rice	PGE13A-50	MODfinder Rapeseed (Canola)	PGE19A-50
MODfinder Potato	PGE20A-50	MODfinder Cotton	PGE21A-50
MODfinder Sugar Beet	PGE22A-50	MODfinder Corn	PGE24A-50
MODfinder Tomato	PGE28A-50	MODfinder Wheat	PGE29A-50
MODfinder Papaya	PGE31A-50		

D. Multiplex* Real-Time PCR kits for botanical impurities detection (50 reactions)

MODfinder MultiENDO 4-plex Corn/Soy/Rapeseed/Cotton	PGE33A-50
MODfinder MultiENDO 2-plex Corn/Soy	PGE43A-50
MODfinder MultiENDO 3-plex Rice/Sugar Beet/Wheat	PGE48A-50
MODfinder MultiENDO 4-plex Flax/Soy/Rapeseed/Cotton	PGE49A-50
MODfinder MultiENDO 2-plex Rapeseed/Cotton	PGE50A-50

E. Multiplex* Real-Time PCR kits for viral particles detection (50 reactions)

MODfinder MultiSCREEN 2-plex CaMV/FMV	PGE39A-50
---------------------------------------	-----------

* 2-plex kits (FAM/HEX); 3-plex kits (FAM/HEX/Cy5); 4-plex (FAM/HEX/Texas Red/Cy5)

Instruction for ordering - MODfinder Identification kits

A. MODfinder singleplex (FAM) Real-Time PCR kits for GM corn event ID (50 reactions)

PGC01A-50	Corn Bt176 (UID SYN-EV176-9)*	PGC02A-50	Corn MON810 (UID MON-00810-6)
PGC03A-50	Corn T25 (UID ACS-ZM003-2)	PGC04A-50	Corn BT11 (UID SYN-BT011-1)
PGC05A-50	Corn NK603 (UID MON-00603-6)	PGC06A-50	Corn CBH-351 (UID ACS-ZM004-3)**
PGC07A-50	Corn GA21 (UID MON-00021-9)	PGC08A-50	Corn DAS1507 (UID DAS-01507-1)
PGC09A-50	Corn MON863 (UID MON-00863-5)*	PGC10A-50	Corn MIR604 (UID SYN-IR604-5)
PGC11A-50	Corn DAS59122 (UID DAS-59122-7)	PGC12A-50	Corn DBT418 (UID DKB-89614-9)**
PGC13A-50	Corn 98140 (UID DP-098140-6)**	PGC14A-50	Corn MIR162 (UID SYN-IR162-4)
PGC16A-50	Corn LY038 (UID REN-00038-3)**	PGC17A-50	Corn Event 3272 (UID SYN-E3272-5)*
PGC19A-50	Corn MON88017 (UID MON-88017-3)	PGC20A-50	Corn MON89034 (UID MON-89034-3)
PGC22A-50	Corn DAS40278 (UID DAS-40278-9)	PGC23A-50	Corn MON87460 (UID MON-87460-4)
PGC24A-50	Corn Event 5307 (UID SYN-05307-1)	PGC25A-50	Corn MON87427 (UID MON-87427-7)
PGC26A-50	Corn VCO-01981-5 (UID VCO-01981-5)**	PGC27A-50	Corn MON87411 (UID MON-87411-9)
PGC29A-50	Corn 4114 (UID DP-004114-3)	PGC30A-50	Corn MZIR098 (UID SYN-00098-3)
PGC31A-50	Corn MZHG0JG (UID SYN-000JG-2)	PGC32A-50	Corn MON87403 (UID MON-87403-1)
PGC42A-50	Corn DP202216 (DP-202216-6)**		

B. 2-plex# MODfinder MultiSCREEN Real-Time PCR kits for GM corn events ID (50 reactions)

PGC28A-50	Corn DAS40278/LY038	PGC33A-50	Corn 4114/5307
PGC34A-50	Corn Bt11/DAS1507	PGC35A-50	Corn GA21/DAS59122
PGC36A-50	Corn MIR162/MIR604	PGC37A-50	MON810/MON87403
PGC38A-50	Corn MON87411/MON87427	PGC39A-50	Corn MON87460/MON88017
PGC40A-50	Corn MON89034/NK603	PGC41A-50	Corn T25/MZHG0JG

C. MODfinder singleplex (FAM) Real-Time PCR kits for GM soy event ID (50 reactions)

PGS01A-50	Soy A2704-12 (UID ACS-GM005-3)	PGS02A-50	Soy GTS 40-3-2 (RoundUp Ready) (UID MON-04032-6)
PGS03A-50	Soy A5547-127 (UID ACS-GM006-4)	PGS04A-50	Soy MON89788 (UID MON-89788-1)
PGS05A-50	Soy DP305423 (UID DP-305423-1)	PGS06A-50	Soy DP356043 (UID DP-356043-5)*
PGS07A-50	Soy MON87701 (UID MON-87701-2)	PGS08A-50	Soy CV127 (UID BPS-CV127-9)
PGS09A-50	Soy MON87705 (UID MON-87705-6)	PGS10A-50	Soy MON87769 (UID MON-87769-7)
PGS11A-50	Soy FG72 (UID MST-FG072-2)	PGS12A-50	Soy DAS44406 (UID DAS-44406-6)
PGS13A-50	Soy DAS68416 (UID DAS-68416-4)	PGS14A-50	Soy DAS81419 (UID DAS-81419-2)
PGS15A-50	Soy MON87708 (UID MON-87708-9)	PGS19A-50	Soy MON87751 (UID MON-87751-7)
PGS20A-50	Soy SYHT0H2 (UID SYN-000H2-5)	PGS33A-50	Soy GMB151 (UID BCS-GM151-6)*

D. Multiplex# MODfinder MultiSCREEN Real-Time PCR kits for GM soy events ID (50 reactions)

PGS17A-50	2-plex Soy CV127/DP305423	PGS22A-50	2-plex Soy SYHT0H2/GMB151
PGS23A-50	2-plex Soy MON87754/MON87751	PGS24A-50	2-plex Soy A2704-12/A5547-127
PGS25A-50	2-plex Soy CV127/MON87769	PGS26A-50	2-plex Soy DAS68416/DAS44406
PGS27A-50	2-plex Soy DP305423/DP356043	PGS28A-50	2-plex Soy GTS 40-3-2/FG72
PGS29A-50	2-plex Soy MON87705/MON89788	PGS30A-50	2-plex Soy MON87708/MON87701
PGS16A-50	3-plex Soy CV127/DP305423/MON87701	PGS18A-50	3-plex Soy A2704/A5547/DP356043
PGS21A-50	3-plex Soy CV127+MON87701/DP305423+DP356043/MON87708+MON87769		
PGS31A-50	4-plex Soy CV127/DP305423/MON87751/MON87708		
PGS32A-50	4-plex Soy MON87701/DP356043/MON87754(fE9)/MON87769		

2-plex kits (FAM/HEX); 3-plex kits (FAM/HEX/Cy5); 4-plex (FAM/HEX/Texas Red/Cy5)

*Events in EU approved with restriction, phasing-out or pending approval – Inquire for availability

**Events non approved in EU – Inquire for availability

E. MODfinder singleplex (FAM) Real-Time PCR kits for GM cotton event ID (50 reactions)

PGT01A-50	Cotton MON531 (UID MON-00531-6)	PGT02A-50	Cotton MON1445 (UID MON-01445-2)
PGT03A-50	Cotton 281-24-236 (UID DAS-24236-5)	PGT04A-50	Cotton 3006-210-23 (UID DAS-21023-5)
PGT05A-50	Cotton LL25 (UID ACS-GH001-3)	PGT06A-50	Cotton MON15985 (UID MON-15985-7)
PGT07A-50	Cotton GHB614 (UID BCS-GH002-5)	PGT08A-50	Cotton MON88913 (UID MON-88913-8)
PGT09A-50	Cotton GHB119 (UID BCS-GH005-8)	PGT10A-50	Cotton T304-40 (UID BCS-GH004-7)
PGT11A-50	Cotton COT-102 (UID SYN-IR102-7)**	PGT12A-50	Cotton MON88701 (UID MON-88701-3)**
PGT13A-50	Cotton DAS81910 (UID DAS-81910-7)**	PGT14A-50	Cotton GHB811 (UID BCS-GH811-4)*

F. MODfinder singleplex (FAM) Real-Time PCR kits for GM rapeseed event ID (50 reactions)

PGZ01A-50	Rapeseed MS8 (UID ACS-BN005-8)	PGZ02A-50	Rapeseed RF3 (UID ACS-BN003-6)
PGZ03A-50	Rapeseed RT73 (UID MON-00073-7)	PGZ04A-50	Rapeseed T45 (UID ACS-BN008-2)
PGZ05A-50	Rapeseed HCN92 (UID ACS-BN007-1)*	PGZ06A-50	Rapeseed RF1 (UID ACS-BN001-4)*
PGZ07A-50	Rapeseed RF2 (UID ACS-BN002-5)*	PGZ08A-50	Rapeseed OXY-235 (UID ACS-BN011-5)**
PGZ09A-50	Rapeseed MS1 (UID ACS-BN004-7)*	PGZ10A-50	Rapeseed DP073496 (UID DP-073496-4)*
PGZ11A-50	Rapeseed MON88302 (UID MON-88302-9)**	PGZ12A-50	Rapeseed MS11 (UID ACS-BN012-7)**

G. MODfinder singleplex (FAM) Real-Time PCR kits for various GM events ID (50 reactions)

PGB01A-50	Sugar beet RURH7-1 (UID KM-000H71-4)	PGF01A-50	Flax FP967 (UID CDC-FL001-2)
PGP01A-50	Potato EH92-527 (UID BPS-25271-9)	PGR01A-50	Rice LL601 (UID BCS-OS003-7)
PGR02A-50	Rice LL62 (UID ACS-OS002-5)	PGR03A-50	Rice BT63 (UID HZU-HH001-9)
PGX01A-50	Salmon AquAdvantage	PGY01A-50	Papaya Huanhong N1
PGY02A-50	Papaya 55-1 (UID CUH-CP551-8)		

MODfinder quantitative kits were developed to quantify the presence of a GMO in a sample according to EU indications

MODfinder portfolio includes purified DNA extracted from traceable reference materials at a defined concentration.

Quantitation of GM content is based on the relative ratio of the copy number of the GM target (detected using the identification PCR reactions) and the copies of the corresponding species detected using the species-specific PCR reactions. Calculation requires two standard curves obtained by using DNA extracts from reference materials.

When testing matrices derived from single ingredients (soy grains, lecithin, corn, corn starch...) a convenient semi-quantitative estimation of GMO contamination can be performed through sample DNA concentration normalization.

A. Quantified Plant DNA GM-free extracts from ISO 17034 validated entities

Each vial contains 120 µL (Conc. 10 ng/ µL) of DNA; inquire for customized formats.

PGE09R-BLK10PPM	Soy	PGE13R-BLK10PPM	Rice	PGE19R-BLK10PPM	Rapeseed
PGE20R-BLK10PPM	Potato	PGE21R-BLK10PPM	Cotton	PGE22R-BLK10PPM	Sugarbeet
PGE24R-BLK10PPM	Corn				

B. Quantified GM Plant DNA extracts from ISO 17034 validated entities

Each vial contains 120 µL (Conc. 10 ng/ µL) of DNA with 1% contamination; inquire for customized formats.

PGB01R	Sugarbeet RURH7-1 (UID KM-000H71-4)	PGS07R	Soy MON87701 (UID MON-87701-2)
PGC01R	Corn Bt176 (UID SYN-EV176-9)	PGS08R	Soy CV127 (UID BPS-CV127-9)
PGC02R	Corn MON810 (UID MON-00810-6)	PGS09R	Soy MON87705 (UID MON-87705-6)
PGC03R	Corn T25 (UID ACS-ZM003-2)	PGS10R	Soy MON87769 (UID MON-87769-7)
PGC04R	Corn BT11 (UID SYN-BT011-1)	PGS11R	Soy FG72 (UID MST-FG072-2)
PGC05R	Corn NK603 (UID MON-00603-6)	PGS12R	Soy DAS44406 (UID DAS-44406-6)
PGC07R	Corn GA21 (UID MON-00021-9)	PGS13R	Soy DAS68416 (UID DAS-68416-4)
PGC08R	Corn DAS1507 (UID DAS-01507-1)	PGS14R	Soy DAS81419 (UID DAS-81419-2)
PGC09R	Corn MON863 (UID MON-00863-5)	PGS15R	Soy MON87708 (UID MON-87708-9)
PGC10R	Corn MIR604 (UID SYN-IR604-5)	PGS19R	Soy MON87751 (UID MON-87751-7)
PGC11R	Corn DAS59122 (UID DAS-59122-7)	PGS20R	Soy SYHT0H2 (UID SYN-000H2-5)
PGC13R	Corn 98140 (UID DP-098140-6)	PGS33R	Soy GMB151 (UID BCS-GM151-6)
PGC14R	Corn MIR162 (UID SYN-IR162-4)	PGT01R	Cotton MON531 (UID MON-00531-6)
PGC17R	Corn Event 3272 (UID SYN-E3272-5)	PGT02R	Cotton MON1445 (UID MON-01445-2)
PGC19R	Corn MON88017 (UID MON-88017-3)	PGT03R	Cotton 281-24-236 (UID DAS-24236-5)
PGC20R	Corn MON89034 (UID MON-89034-3)	PGT04R	Cotton 3006-210-23 (UID DAS-21023-5)
PGC22R	Corn DAS40278 (UID DAS-40278-9)	PGT05R	Cotton LL25 (UID ACS-GH001-3)
PGC23R	Corn MON87460 (UID MON-87460-4)	PGT06R	Cotton MON15985 (UID MON-15985-7)
PGC24R	Corn Event 5307 (UID SYN-05307-1)	PGT07R	Cotton GHB614 (UID BCS-GH002-5)
PGC25R	Corn MON87427 (UID MON-87427-7)	PGT08R	Cotton MON88913 (UID MON-88913-8)
PGC26R	Corn VCO-01981-5 (UID VCO-01981-5)	PGT09R	Cotton GHB119 (UID BCS-GH005-8)
PGC27R	Corn MON87411 (UID MON-87411-9)	PGT10R	Cotton T304-40 (UID BCS-GH004-7)
PGC29R	Corn 4114 (UID DP-004114-3)	PGT11R	Cotton COT-102 (UID SYN-IR102-7)
PGC30R	Corn MZIR098 (UID SYN-00098-3)	PGT12R	Cotton MON88701 (UID MON-88701-3)
PGC31R	Corn MZHG0JG (UID SYN-000JG-2)	PGT13R	Cotton DAS81910 (UID DAS-81910-7)
PGC32R	Corn MON87403 (UID MON-87403-1)	PGT14R	Cotton GHB811 (UID BCS-GH811-4)
PGC42R	Corn DP202216 (UID DP-202216-6)	PGZ01R	Rapeseed MS8 (UID ACS-BN005-8)
PGP01R	Potato EH92-527 (UID BPS-25271-9)	PGZ02R	Rapeseed RF3 (UID ACS-BN003-6)
PGP02R	Potato AM04-1020 (UID BPS-A1020-5)	PGZ03R	Rapeseed GT73 (UID MON-00073-7)
PGP03R	Potato AV43-6-G7 (UID AVE-436G7-1)	PGZ04R	Rapeseed T45 (UID ACS-BN008-2)
PGP04R	Potato PH05-026-0048 (UID BPS-PH048-1)	PGZ05R	Rapeseed TOPAS 19 2 (UID ACS-BN007-1)
PGR02R	Rice LL62 (UID ACS-OS002-5)	PGZ06R	Rapeseed RF1 (UID ACS-BN001-4)
PGS01R	Soy A2704-12 (UID ACS-GM005-3)	PGZ07R	Rapeseed RF2 (UID ACS-BN002-5)
PGS02R	Soy GTS 40-3-2 (UID MON-04032-6)	PGZ09R	Rapeseed MS1 (UID ACS-BN004-7)
PGS03R	Soy A5547-127 (UID ACS-GM006-4)	PGZ10R	Rapeseed DP073496 (UID DP-073496-4)
PGS04R	Soy MON89788 (UID MON-89788-1)	PGZ11R	Rapeseed MON88302 (UID MON-88302-9)
PGS05R	Soy DP305423 (UID DP-305423-1)	PGZ12R	Rapeseed MS11 (UID ACS-BN012-7)
PGS06R	Soy DP356043 (UID DP-356043-5)	PGZ13R	Rapeseed LBFLFK (UID BPS-BFLFK-2)

Each kit contains: event specific primers and probe mix (50 tests, 30 µl reaction end volume); taxon specific primers and probe mix (50 tests, 30 µl reactions end volume); standard solutions at scalar target GM concentration; taxon DNA standard solutions at scalar concentration; negative control. Standard solutions are obtained extracting lot traceable reference materials using Ion-Force DNA Extraction kit.

A. MODfinder Real-Time PCR kits for GM corn quantification (50 reactions)

PGC01Q-50	Corn Bt176 (UID SYN-EV176-9)*	PGC02Q-50	Corn MON810 (UID MON-00810-6)
PGC03Q-50	Corn T25 (UID ACS-ZM003-2)	PGC04Q-50	Corn BT11 (UID SYN-BT011-1)
PGC05Q-50	Corn NK603 (UID MON-00603-6)	PGC07Q-50	Corn GA21 (UID MON-00021-9)
PGC08Q-50	Corn DAS1507 (UID DAS-01507-1)	PGC09Q-50	Corn MON863 (UID MON-00863-5)*
PGC10Q-50	Corn MIR604 (UID SYN-IR604-5)	PGC11Q-50	Corn DAS59122 (UID DAS-59122-7)
PGC13Q-50	Corn 98140 (UID DP-098140-6)**	PGC14Q-50	Corn MIR162 (UID SYN-IR162-4)
PGC17Q-50	Corn Event 3272 (SYN-E3272-5)*	PGC19Q-50	Corn MON88017 (UID MON-88017-3)
PGC20Q-50	Corn MON89034 (UID MON-89034-3)	PGC22Q-50	Corn DAS40278 (UID DAS-40278-9)
PGC23Q-50	Corn MON87460 (UID MON-87460-4)	PGC24Q-50	Corn Event 5307 (UID SYN-05307-1)
PGC25Q-50	Corn MON87427 (UID MON-87427-7)	PGC26Q-50	Corn VCO-01981 (VCO-01981-5)**
PGC27Q-50	Corn MON87411 (UID MON-87411-9)	PGC29Q-50	Corn 4114 (UID DP-004114-3)
PGC30Q-50	Corn MZIR098 (UID SYN-00098-3)	PGC31Q-50	Corn MZHGOJG (UID SYN-000JG-2)
PGC32Q-50	Corn MON87403 (UID MON-87403-1)	PGC42Q-50	Corn DP202216 (DP-202216-6)**

B. MODfinder Real-Time PCR kits for GM soy quantification (50 reactions)

PGS01Q-50	Soy A2704-12 (UID ACS-GM005-3)	PGS02Q-50	Soy GTS 40-3-2 (UID MON-04032-6)
PGS03Q-50	Soy A5547-127 (UID ACS-GM006-4)	PGS04Q-50	Soy MON89788 (UID MON-89788-1)
PGS05Q-50	Soy DP305423 (UID DP-305423-1)	PGS06Q-50	Soy DP356043 (UID DP-356043-5)*
PGS07Q-50	Soy MON87701 (UID MON-87701-2)	PGS08Q-50	Soy CV127 (UID BPS-CV127-9)
PGS09Q-50	Soy MON87705 (UID MON-87705-6)	PGS10Q-50	Soy MON87769 (UID MON-87769-7)
PGS11Q-50	Soy FG72 (UID MST-FG072-2)	PGS12Q-50	Soy DAS44406 (UID DAS-44406-6)
PGS13Q-50	Soy DAS68416 (UID DAS-68416-4)	PGS14Q-50	Soy DAS81419 (UID DAS-81419-2)
PGS15Q-50	Soy MON87708 (UID MON-87708-9)	PGS19Q-50	Soy MON87751 (UID MON-87751-7)
PGS20Q-50	Soy SYHT0H2 (UID SYN-000H2-5)	PGS33Q-50	Soy GMB151 (UID BCS-GM151-6)*

C. MODfinder Real-Time PCR kits for GM cotton quantification (50 reactions)

PGT01Q-50	Cotton MON531 (UID MON-00531-6)	PGT02Q-50	Cotton MON1445 (UID MON-01445-2)
PGT03Q-50	Cotton 281-24-236 (UID DAS-24236-5)	PGT04Q-50	Cotton 3006-210-23 (UID DAS-21023-5)
PGT05Q-50	Cotton LL25 (UID ACS-GH001-3)	PGT06Q-50	Cotton MON15985 (UID MON-15985-7)
PGT07Q-50	Cotton GHB14 (UID BCS-GH002-5)	PGT08Q-50	Cotton MON88913 (UID MON-88913-8)
PGT09Q-50	Cotton GHB119 (UID BCS-GH005-8)	PGT10Q-50	Cotton T304-40 (UID BCS-GH004-7)
PGT11Q-50	Cotton COT-102 (UID SYN-IR102-7)**	PGT12Q-50	Cotton MON88701 (UID MON-88701-3)**
PGT13Q-50	Cotton DAS81910 (UID DAS-81910-7)**	PGT14Q-50	Cotton GHB811 (UID BCS-GH811-4)*

D. MODfinder Real-Time PCR kits for GM rapeseed quantification (50 reactions)

PGZ01Q-50	Rapeseed MS8 (UID ACS-BN005-8)	PGZ02Q-50	Rapeseed RF3 (UID ACS-BN003-6)
PGZ03Q-50	Rapeseed GT73 (UID MON-00073-7)	PGZ04Q-50	Rapeseed T45 (UID ACS-BN008-2)
PGZ05Q-50	Rapeseed Topas19/2 (UID ACS-BN007-1)*	PGZ06Q-50	Rapeseed RF1 (UID ACS-BN001-4)*
PGZ07Q-50	Rapeseed RF2 (UID ACS-BN002-5)*	PGZ09Q-50	Rapeseed MS1 (UID ACS-BN004-7)*
PGZ10Q-50	Rapeseed DP073496 (UID DP-073496-4)*	PGZ11Q-50	Rapeseed MON88302 (UID MON-88302-9)
PGZ12Q-50	Rapeseed MS11 (UID ACS-BN012-7)**		

E. MODfinder Real-Time PCR kits for various GM events GM quantification (50 reactions)

PGB01Q-50	Sugarbeet H7-1 (KM-000H71-4)
-----------	------------------------------

*Events in EU approved with restriction, phasing-out or pending approval – Inquire for availability

**Events non approved in EU – Inquire for availability



INGREDIENTS AUTHENTICITY | CHEMICAL RESIDUES
MICROBIOLOGY | MYCOTOXINS
ALLERGENS | GMO

Generon S.p.A.

Via San Geminiano 4 – San Prospero (MO) – Italy - ☎ +39 059 8637161 📞 +39 059 7353024

✉ marketing@generon.it 🌐 www.generon.it